



DK50 D DK50 DM *event*

МЕДИЦИНСКИЕ КОМПРЕССОРЫ

Для подачи чистого сжатого отфильтрованного воздуха высокого качества.

Качество воздуха соответствует стандарту EN DIN 12021.

Индикация высокой температуры внутри шкафа компрессора.

Встроенный сепаратор воды понижает точку росы ниже температуры окружающей среды.

Безмасляный компрессор с термозащитой. Реле давления автоматически включает компрессор.

Давление сжатого воздуха на выходе регулируется фильтром регулятором и отображается манометром.

- **ПЕРЕДВИЖНЫЙ**
- **НАДЕЖНЫЙ**
- **ТИХИЙ**



Тех. параметры	DK50 D	DK50 DM	
Номинальное напряжение/ частота (*)	230/50, 230/60 120/60 100/50..60	230/50, 230/60 120/60	В/Гц
Производительность компрессора	40	60	Л.мин ⁻¹ / 3 бар
Пик-флюу при давлении 3 бар		200	Л.мин ⁻¹
Объем рисивера	5		л
Степень фильтрации	5		микронов
Понижение точки росы при 3 бар/ 40 Lit.min ⁻¹ , 20°C	5°C градусов ниже температуры окружающей среды	-	
Понижение точки росы при 3 бар/ 60 Lit.min ⁻¹ , 20°C	-	5°C градусов ниже температуры окружающей среды	
Коннекторы на выходе	RECTUS 21KS KS08 MPC		
Режим работы	Бесперерывный - C1		
Уровень шума	≤ 49	≤ 50	дБ [A]
Размеры (ш x л x в)	400 x 430 x 780	400 x 500 x 780	мм
Весъ	46	59	кг

* под заказ - требуемое напряжение просим указать в заявке



DK50 D

DK50 DM

event

MEDICAL COMPRESSORS

For the delivery of dry medical grade quality compressed air.

Air quality acc.to EN DIN 12021.

Indication of high temperature inside the compressor's box.



Integrated water separator keeps the dew point below ambient temperature.

Oil free compressors with temperature protection.

Automatic pressure switching of the compressor to a small air tank.

Output air pressure set by filter regulator and indicated by manometer.

- PROMPTLY
- RELIABLE
- QUIET

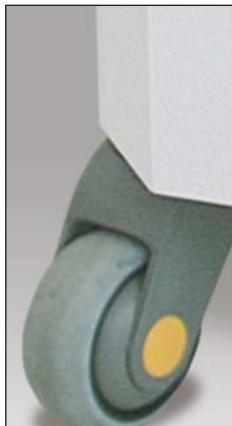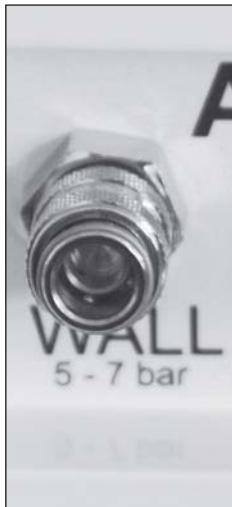

Technical specification	DK50 D	DK50 DM	
Nominal voltage and frequency (*)	230/50, 230/60 120/60 100/50..60	230/50, 230/60 120/60	V/Hz
Compressor efficiency	40	60	$\text{l} \cdot \text{min}^{-1} / 3 \text{ bar}$
Peak flow at pressure 3 bar		200	$\text{l} \cdot \text{min}^{-1}$
Air tank capacity	5		lit.
Filtration of air to	5		μm
Dew point depression at 3 bar/ 40 $\text{lit} \cdot \text{min}^{-1}$, 20°C	5°C under the ambient temperature	-	
Dew point depression at 3 bar/ 60 $\text{lit} \cdot \text{min}^{-1}$, 20°C	-	5°C under the ambient temperature	
Outlet connection	RECTUS 21KS KS08 MPC		
Mode of the operation	Continuous - S1		
Noise level	≤ 49	≤ 50	dB [A]
Dimensions (w x d x h)	400 x 430 x 780	400 x 500 x 780	mm
Weight	46	59	kg

* please state voltage and frequency in the order